
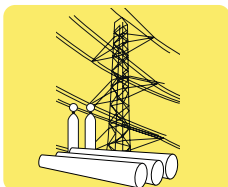


# Emergency User Guide

The SAM 360P will only administer a shock if it is needed and will **NOT** harm a person.  
Turn on  the SAM 360P and follow the voice prompts.



## 1. Remove danger

Remove any source of danger or move the victim away from danger. Ensure your own safety.



**Caution: You must use the SAM 360P at least 2 m (6 feet) away from all radio frequency devices. Alternatively switch off the equipment causing the electromagnetic interference.**



## 2. Check for a response

Shake the victim by the shoulders. Speak loudly to the victim.



## 3. Check the airway

Check for breathing. Open the airway with a head tilt-chin lift if necessary.

Continued overleaf...

#### 4. Get help


Call the emergency services. Get other people to help you. Get the SAM 360P.

You do not need to remove the SAM 360P from its soft carry case to use it.



#### 5. Perform CPR

Push hard and fast at a rate of at least 100 beats per minute, 5 to 6 cm deep, on victim's chest. 30 compressions followed by 2 breaths.

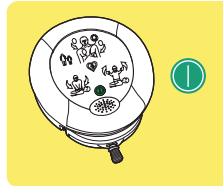


#### 6. Turn on the SAM 360P

Turn on the SAM 360P.

Depending on the Pad-Pak currently installed, the SAM 360P will give a voice prompt 'adult patient' or 'child patient'. If necessary, change the Pad-Pak to the correct type.




**Warning: If a Pediatric-Pak, or an alternative suitable defibrillator is not available, you may use an adult Pad-Pak to treat a child.**

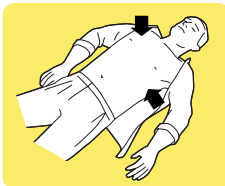


## 7. Bare the chest area

Remove the victim's clothes from the chest area. Remove any metal (bras and jewellery) from the chest area. Make sure the victim's chest is dry. If the victim has a lot of chest hair, shave the part of the chest where you will place the electrodes.

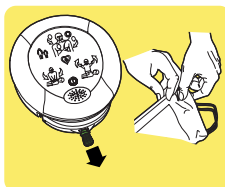


## 8. Open the electrode pouches

Pull on the green tab to open the electrode drawer and tear open the electrode pouch.

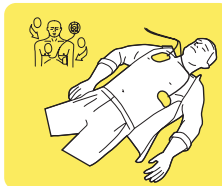

## 9. Place the electrode pads

Follow the voice prompts. Peel the backing off the electrode pads and place them on the chest as shown.

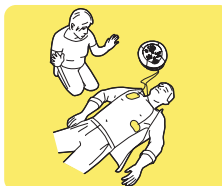

## 10. Do NOT touch the victim

Make sure no one touches the victim while the SAM 360P is analysing or delivering a shock.

The SAM 360P will automatically administer a shock if it is needed. A voice prompt will tell you what to do.

Continued

overleaf...

## 11. Continue

Continue following the voice prompts. The SAM 360P will tell you to continue with CPR if it is necessary.

The metronome will help you to maintain correct CPR. The SAM 360P will continue to give shocks if necessary, but not without warning you first.

If the victim recovers, leave the electrode pads connected and continue to follow the voice prompts until the emergency services arrive.



### Worldwide Headquarters:

HeartSine Technologies Inc. / 121 Friends Lane / Suite 400 / Newtown, PA 18940 / USA  
Tel: 1.215.860.8100 / Fax: 1.215.860.8192

### Manufactured by:

HeartSine Technologies Ltd. / Canberra House / 203 Airport Road West / Belfast BT3 9ED / Northern Ireland  
Tel: +44 (0)28 9093 9400 / Fax: +44 (0)28 9093 9401

[www.heartsine.com](http://www.heartsine.com)  
[support@heartsine.com](mailto:support@heartsine.com)

H037-019-226-0 February 2014 (UK English)

